

### Automated import of patient clinical and risk factors.

PenRad facilitates a one-step process to import as a constructed phrase into the clinical section of the report; Patient clinical and risk factors, family breast cancer history, plus BRCA and Tyrer Cuzick values.

Clinical can also include macros, free text, previous clinical statements, and statements received with order.

To access the clinical section of report tap the **Clinical Text** button. Button is available on Preliminary exam screen, PreExam screen, each of the Examination main screens, and on approve screens.

To import clinical macros, select section from Phrase Category by tapping on category, then tap on desired phrases to incorporate (multiple macros can be used in conjunction). See below sample.

To import patient risk factors, family breast cancer history, BRCA and Tyrer Cuzick generated phrase, tap **Insert Risk History** button. See below sample.



To edit clinical phrase, edit in **Clinical Text** window. Macro and imports insert at insert cursor.

To view in patient's risk history in persistent window during exam, tap **Risk History Summary** button.

To view, use patient's previous incorporated clinical statement, or clinical provided with order, tap **Import Clinical** button. To reuse statement (can be edited), tap to select and tap **Select** button.

Sample below incorporates macro and patient risk and history summary.

Tap on macro, macro text inserted, then tap Risk History Summary button, summary text inserted.

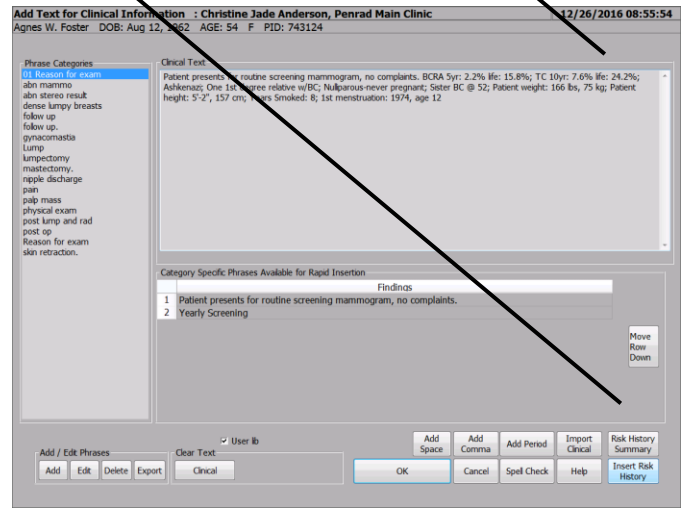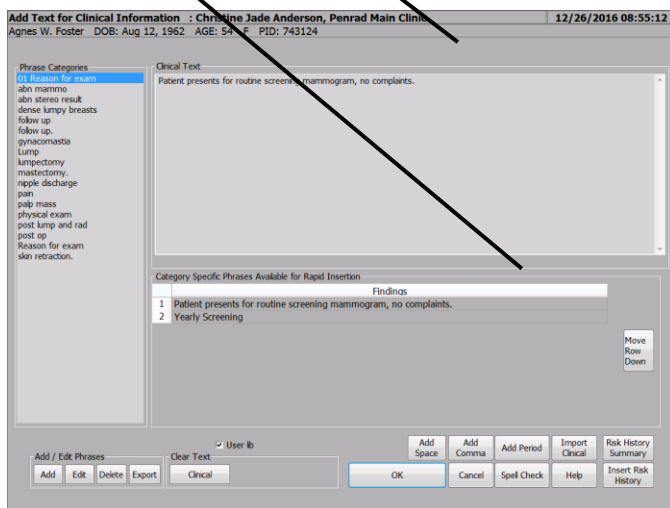

Interested in one-on-one training or a refresher? This insures your facility is incorporating all the latest technology, thus maximizing efficiency and economics. Give us a call at 763-475-3388.

**Recommendation for product development?**

sales@penrad.com | 763.475.3388. Thank you for your business.

© 2017 PenRad Technologies, Inc. 730-23-00401 AA